


Tuff

With their three thread zones, Tuff implants have been uniquely designed according to the anatomy of the bone structure. The lower V-shape thread zone enables self-tapping. The middle zone square type thread is used for compressing cancellous bone and helping achieve maximum BIC. The micro thread on the upper zone adds stability and prevents crestal bone loss.

Material: Titanium (Ti6Al4V ELI)


Treatment: RBM

Available in two neck textures (See page 8 for more information)

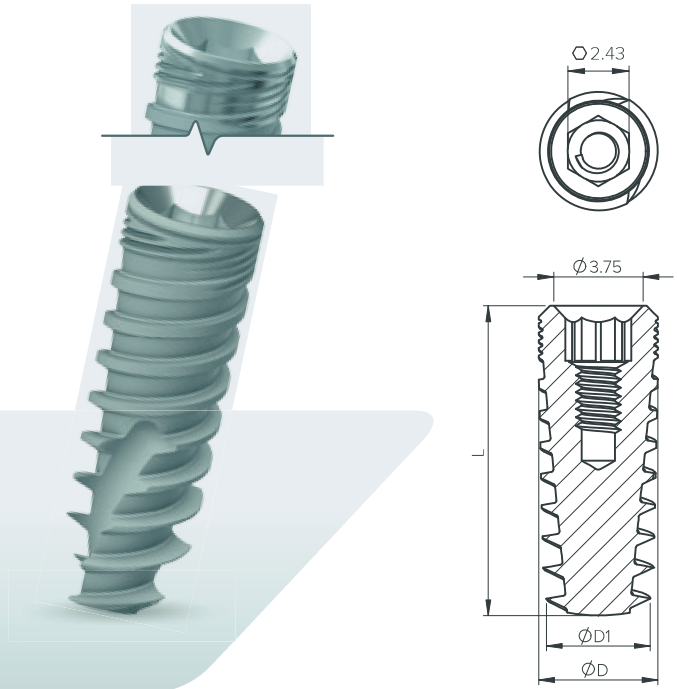


Drill Sequence*

Bone type D1, D2

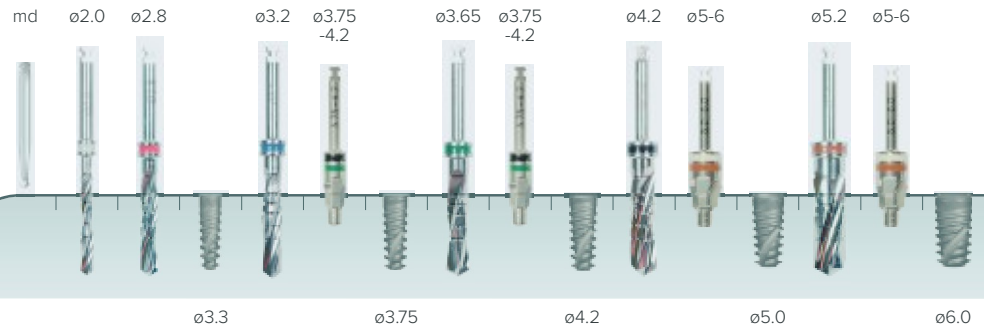





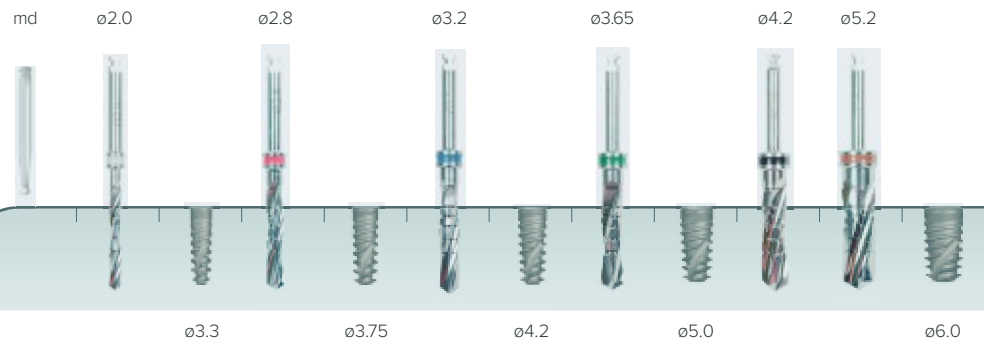

Drill Sequence*

Bone type D3, D4



* The proposed procedure is only a recommendation and should not replace the doctor's judgment. The implants may be placed in immediate function when good primary stability has been achieved and with appropriate occlusal loading.

| Implants | Ø D (mm) | Ø D1 (mm) | L (mm) | RBM Treated Neck | Machined Neck |
|---|----------|-----------|----------|------------------|---------------|
|  | 3.3 | 2.6 | 8 | NM-F3308 | NMSF3308 |
| | | | 10 | NM-F3310 | NMSF3310 |
| | | | 11.5 | NM-F3311 | NMSF3311 |
| | | | 13 | NM-F3313 | NMSF3313 |
| | | | 16 | NM-F3316 | NMSF3316 |
|  | 3.75 | 3.1 | 6 | NM-F3706 | NMSF3706 |
| | | | 8 | NM-F3708 | NMSF3708 |
| | | | 10 | NM-F3710 | NMSF3710 |
| | | | 11.5 | NM-F3711 | NMSF3711 |
| | | | 13 | NM-F3713 | NMSF3713 |
| | | | 16 | NM-F3716 | NMSF3716 |
| | | | 18 | NM-F3718 | NMSF3718 |
|  | 4.2 | 3.5 | 6 | NM-F4206 | NMSF4206 |
| | | | 8 | NM-F4208 | NMSF4208 |
| | | | 10 | NM-F4210 | NMSF4210 |
| | | | 11.5 | NM-F4211 | NMSF4211 |
| | | | 13 | NM-F4213 | NMSF4213 |
| | | | 16 | NM-F4216 | NMSF4216 |
| | | | 18 | NM-F4218 | NMSF4218 |
| | | | 20 | NM-F4220 | NMSF4220 |
| | | | 22 | NM-F4222 | NMSF4222 |
|  | 5.0 | 4.5 | 6 | NM-F5006 | NMSF5006 |
| | | | 8 | NM-F5008 | NMSF5008 |
| | | | 10 | NM-F5010 | NMSF5010 |
| | | | 11.5 | NM-F5011 | NMSF5011 |
| | | | 13 | NM-F5013 | NMSF5013 |
| | | | 16 | NM-F5016 | NMSF5016 |
|  | 6.0 | 5.2 | 6 | NM-F6006 | NMSF6006 |
| | | | 8 | NM-F6008 | NMSF6008 |
| | | | 10 | NM-F6010 | NMSF6010 |
| | | | 11.5 | NM-F6011 | NMSF6011 |
| | | | 13 | NM-F6013 | NMSF6013 |
| | | | 16 | NM-F3308 | NMSF3308 |
| Cover Screw Included with all implants  | | | NM-S5023 | | |